


# Corticosteroids by LC-MS/MS

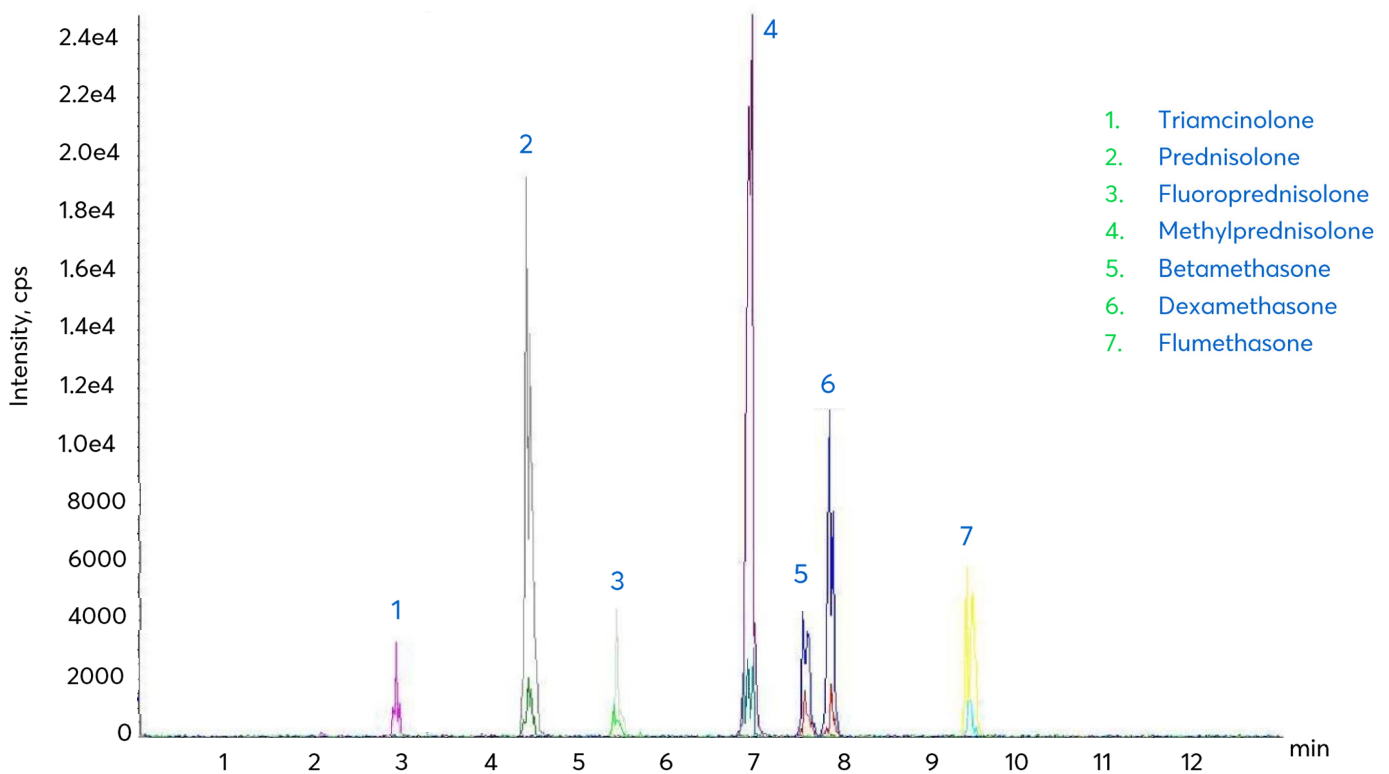

# Method Details

## CONDITIONS

Column: ACE 3 C18-PFP  
 Dimensions: 150 x 2.1 mm  
 Part Number: ACE-1110-1502  
 Mobile Phase: A: 0.1% formic acid in H<sub>2</sub>O  
 B: MeCN/0.1% formic acid in H<sub>2</sub>O (90:10 v/v)

Gradient:

Time (mins)	%B
0.0	30
14.0	50
17.0	95
20.0	30

Flow Rate: 0.3 mL/min  
 Injection: 25 µL  
 Temperature: 15 °C  
 Detection: Turbospray, MRM

## ORDERING TABLE

Product	Details	Size	Part number
Avantor® ACE® 3 C18-PFP	HPLC Column	150 x 2.1 mm	ACE-1110-1502

Reproduced with permission of Laboratorio Salud Publica de Burgos, Spain