


# Application note #2290

## Phenol and Phenoxy Acid Herbicides

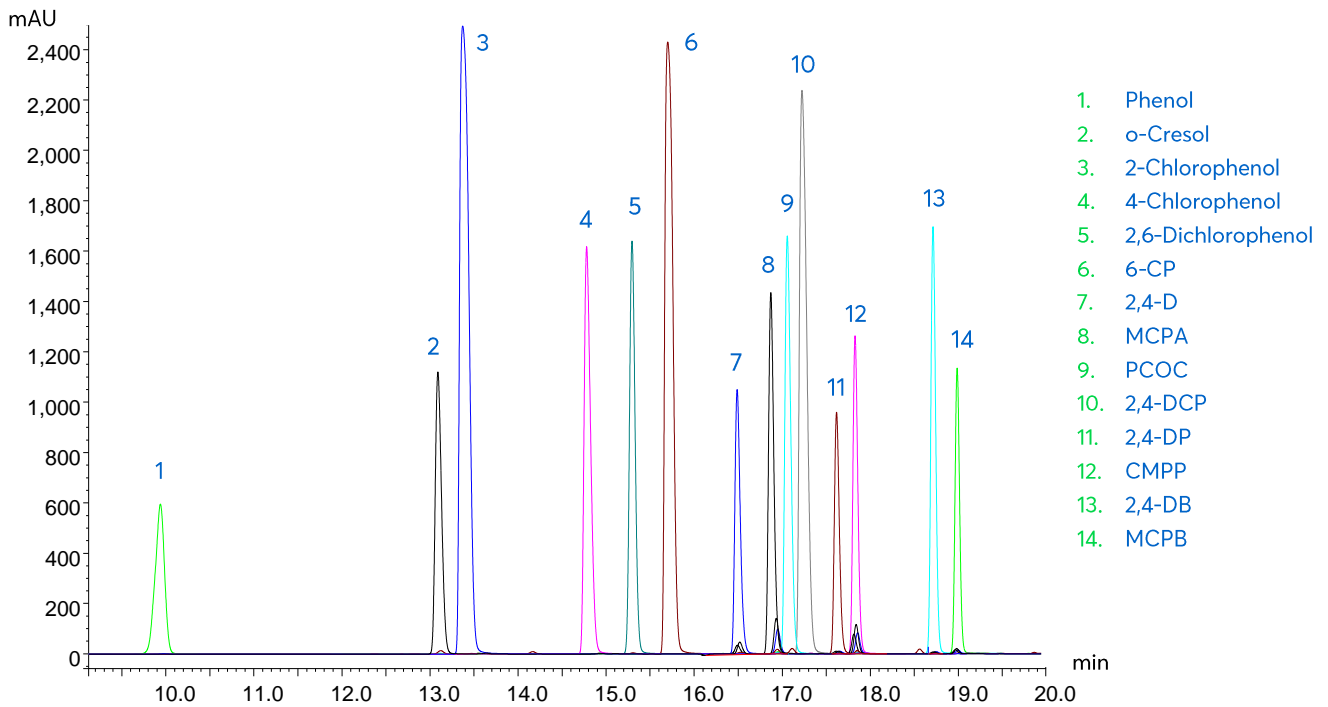

# Method Details

## CONDITIONS

Column: Avantor® ACE® 3 C18-PFP  
 Dimensions: 150 x 4.6 mm  
 Mobile Phase: A: 0.1% formic acid in H<sub>2</sub>O  
 B: MeOH

Time (mins)	% B
0.0	10
20.0	100

Flow Rate: 1 mL/min  
 Injection: 10 µL  
 Temperature: 35 °C  
 Detection: UV, 280 nm

## ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® 3 C18-PFP	HPLC Column	150 x 4.6 mm	ACE-1110-1546

Reproduced with permission of Nufarm UK Ltd